




Features to include in GUI:

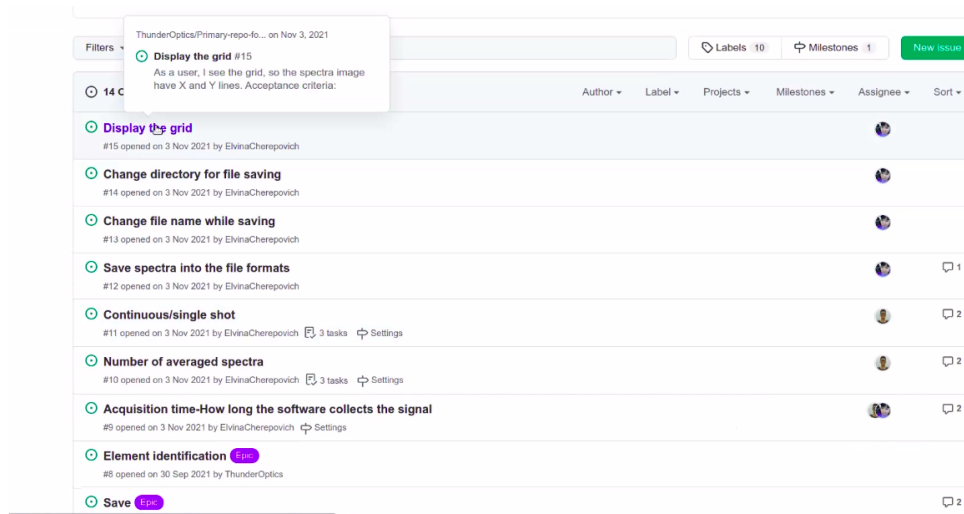

Display the grid #15

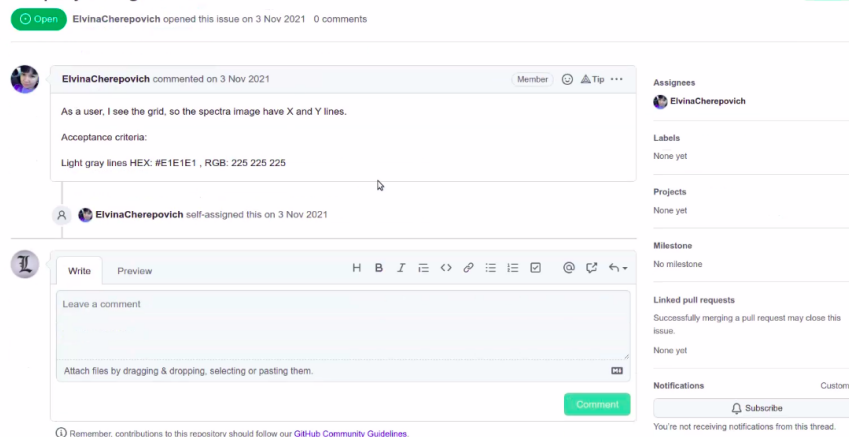

Change file name while saving #13

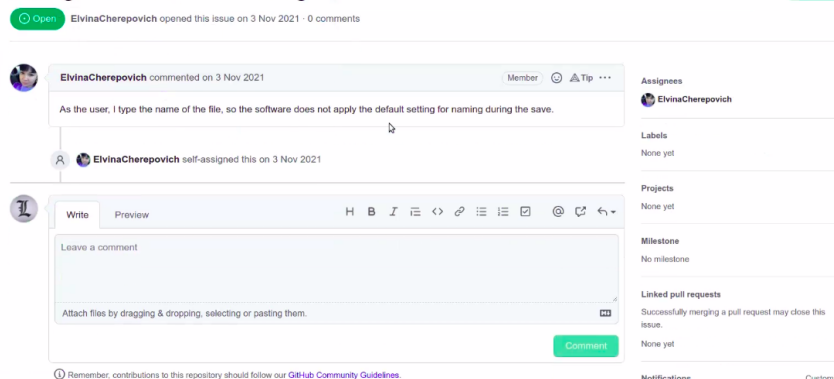

Save spectra into the file formats #12

[Open](#) ElvinaCherepovich opened this issue on 3 Nov 2021 · 1 comment



ElvinaCherepovich commented on 3 Nov 2021

Member



As the user, I choose the format to save, so the software saves spectra into the file of a specific format.

Acceptance criteria:

File formats: excel, dat, pdf.

The default folder is Documents\spectra

The default name of the file is DDMYYYY-N+1. N - is incremental.



ElvinaCherepovich self-assigned this on 3 Nov 2021



ElvinaCherepovich mentioned this issue on 3 Nov 2021

[Change directory for file saving #14](#)

[Open](#)



pm1606 commented on 10 Nov 2021

Member



In progress

Assignees

ElvinaCherepovich

Labels

None yet

Projects

None yet

Milestone

No milestone

Linked pull requests

Successfully merging a pull request may close this issue.

None yet

Notifications

Customize

[Subscribe](#)

Continuous/single shot #11

[Open](#) [3 tasks](#) ElvinaCherepovich opened this issue on 3 Nov 2021 · 2 comments



ElvinaCherepovich commented on 3 Nov 2021

Member



This parameter includes two user stories:

1. As the user, I set the parameter "Continuous" shot, so the software takes spectra nonstop with default intervals.
2. As the user, I set the parameter "Single" shot, so the software takes spectra one time.

Acceptance Criteria:

- ☐ The default value is a Single shot.
- ☐ Software delete a previous shot.
- ☐ The interval for continuous mode is 250ms and cannot be changed.

Consideration:

- Continuous mode allows seeing changes of spectra in real-time if the user changes the intensity of a light source, plays with other settings like acquisition time, and so on.



ElvinaCherepovich added this to the **Settings** milestone on 3 Nov 2021



ElvinaCherepovich assigned pm1606 on 3 Nov 2021



pm1606 commented on 5 Nov 2021

Member



In progress

Assignees

pm1606

Labels

None yet

Projects

None yet

Milestone

Settings

Linked pull requests

Successfully merging a pull request may close this issue.

None yet

Notifications

Customize

[Subscribe](#)

You're not receiving notifications from this thread.

2 participants



Acquisition time-How long the software collects the signal #9

[Open](#) ElvinaCherepovich opened this issue on 3 Nov 2021 · 2 comments



ElvinaCherepovich commented on 3 Nov 2021 · edited

Member



As the user, I set the acquisition time, so during this time, the spectrometer acquires a signal (collect light)

Acceptance Criteria:

1. Default parameter is 100ms.
2. Each step increases or decreases the parameter to 100ms.
3. Unit is ms (milliseconds) until it reaches 900ms.
4. Unit is s (seconds) if the value is more than 900ms s (seconds).
5. Parameters are from 100ms to 10s.
6. Parameters can be scrolled with the mouse.

Considerations:

1. This parameter defines how intense will be the detected signal.
2. 1 s = 1000ms.
3. Unit is ms (milliseconds) until it reaches 1000ms.
4. Unit is s (seconds) if the value is more than 900ms s (seconds).
- 5.



ElvinaCherepovich added this to the **Settings** milestone on 3 Nov 2021

Assignees

ElvinaCherepovich

pm1606

Labels

None yet

Projects

None yet

Milestone

Settings

Linked pull requests

Successfully merging a pull request may close this issue.

None yet

Notifications

Customize

[Subscribe](#)

Number of averaged spectra #10

Open

3 tasks

ElvinaCherepovich opened this issue on 3 Nov 2021 · 2 comments

How does it work?



ElvinaCherepovich commented on 3 Nov 2021

Member



As the user, I set the number of average spectra, so the spectrometer acquires a number of spectra and makes an overage.

Acceptance Criteria:

- ☐ Default is 1.
- ☐ Optional values are from 1 to 100.
- ☐ User type or scroll the parameter.

Consideration:

1. This parameter helps reduce the noise.
2. The number of spectra is also an important parameter of the acquisition.

Spectrum is a set of the values of intensity of each pixel in the linear detector. So, it is a set of numbers (163, 89765, 87, ...). Each number is the intensity of the corresponding pixel. The intensity is the current that each pixel sends, it corresponds also to the amount of light the pixel receives.



ElvinaCherepovich added this to the **Settings** milestone on 3 Nov 2021



ElvinaCherepovich assigned pm1606 on 3 Nov 2021



pm1606 commented on 5 Nov 2021

Member



Assignees

pm1606

Labels

None yet

Projects

None yet

Milestone

Settings

Linked pull requests

Successfully merging a pull request may close this issue.

None yet

Notifications

Customize

Subscribe

You're not receiving notifications from this thread.

2 participants